

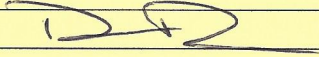
**EQUIPMENT DETAILS:**

Equipment Type/Model:	HEITKA ROTINA 380R
Catalogue Code:	1706
Serial Number:	000 2188 - 04
Warranty Period:	2 years

**Checks Completed**

	Pass
Instrument operational	✓
Test with appropriate rotor where available	✓
Electrically tested (Earth <1 Ω, Insulation >1 MΩ in resistance)	✓
Correct plug fitted, insulated pins	✓
Exterior condition – clean, tidy, no damage	✓
Interior Condition – clean, tidy, no damage	✓
All manuals and accessories included	✓
LGA Equipment Sticker applied	✓
P1 – 3000 g : 10 min : 20 deg : Set Rotor 174 mm for 1754 rotor	✓
P2 – 3000 g : 10 min : 4 deg : Set Rotor 174 mm for 1754 rotor	
P3 – Other rotor as required	
P4 – Other Rotor as required	
Put in Lock 1 Mode	✓
Transport Device : Remove for remote delivery	✓

Parameter	Set Point	Display	Actual	Pass
Speed	3000	3000	2997	✓

Tested by:	Drew Dixon	Signature:	
Date:	27.6.14.		

**\*\*Please note: all test equipment used for this check are traceable to NATA calibrated instruments.**

Test Equipment					
Parameters	Meter	Serial #	Certificate #	Valid To	Test Lab Cert #
Electronic	Fluke 179	21530350	28753-1/1	23/11/2013	15257
Electrical	Kyoritsu 6201A	W0316897	28753-2/1	23/11/2013	15257
Speed	Extech 461920	12079763	28821-1/1	29/11/2013	15257
Temperature	Fluke 179	21530350	28753-1/1	23/11/2013	15257
Timer	iphone based timer synchronised with central atomic time clock server				ptbtime1.ptb.de